


How to Easily Name and Store DSA800 Screen-Shot Files

DSA800 firmware version 00.01.12 and later has the ability to save a prefix for use with names of screen-shot file saves. This can make entering a file name for saved screen-shots much easier.

I. Add a Screen-Shot File Name Prefix:


1. Press STORE > Down arrow to page 2 > Select Name Prefix

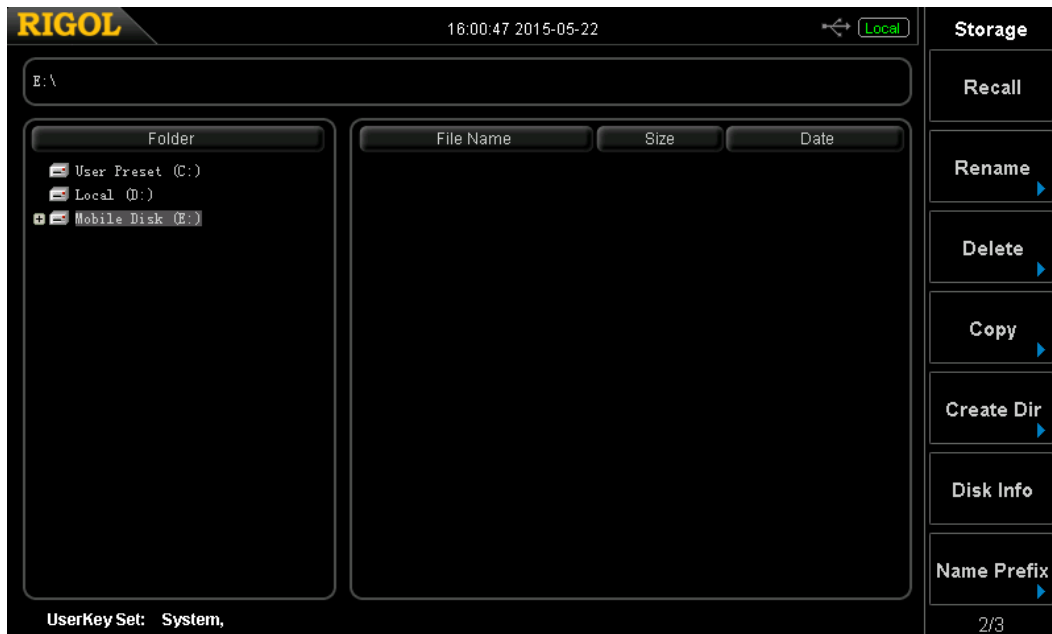

2. Select Edit Prefix.



3. Use the keypad to enter the prefix name. Use the name **TEST** (followed with one blank space) for this example. *Note: The 'file name' shown below doesn't have a blank space after TEST.*

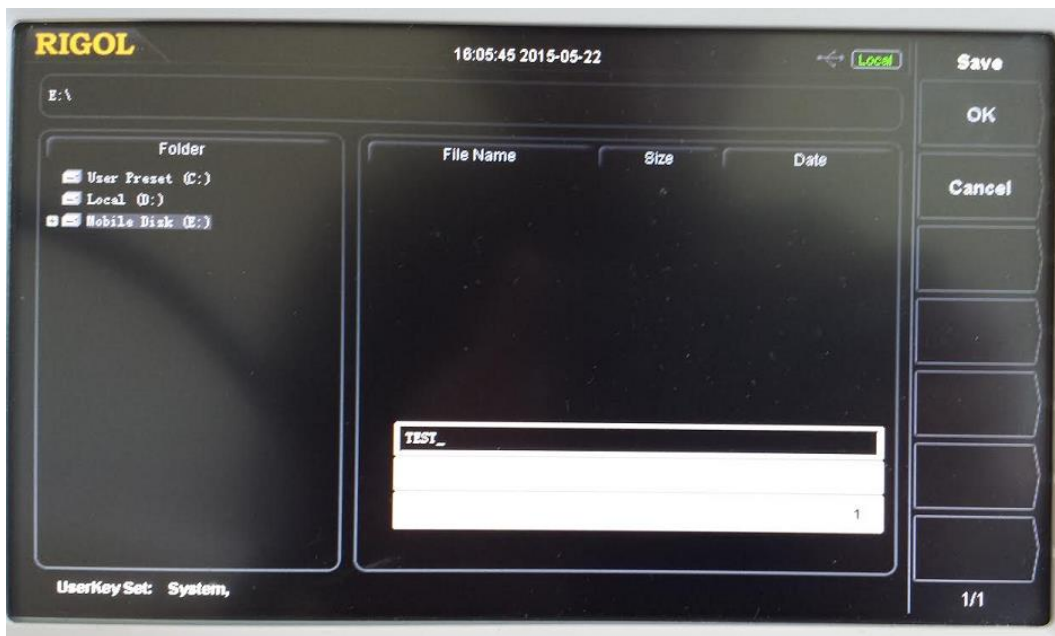

4. **Note:** You can leave the text box in the numeric control mode (with '1' shown in lower corner) when you close the text window. It will appear in the same state the next time you save a file. You can then simply press 1.2.3.. etc. to increment the file name's suffix.
5. Press OK to save the file name prefix TEST (with one blank space behind it).
6. Press Cancel to return to the Name Prefix page.
7. Press Prefix Edit to select ON. This makes the prefix name TEST active for your screen-shot saves. *Note: Prefix Edit can be temporarily turned OFF anytime for names without the TEXT prefix.*

II. Using the File Name Prefix **TEST** With Saved Screen-Shots:

1. Select a 'screen display' to save as a screen-shot on a USB Flash Drive.
2. Insert a USB Flash Drive (FAT32 format recommended).
3. Press the Printer icon button (center, left side of LCD) to save the screen to the USB drive.
4. Enter the suffix to the prefix name TEST. i.e. Add 12 for the name 'TEST 12'